




**STORZ & BICKEL**




***Volcano***<sup>®</sup>  
**CLASSIC**

**INSTRUCTIONS FOR USE (CONCISE)**

## FUNCTIONAL ELEMENTS

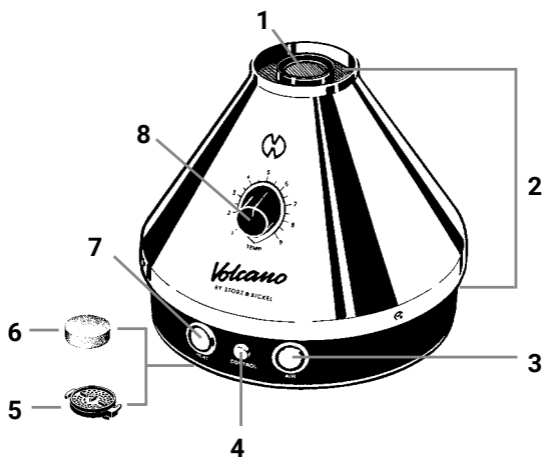

-  **Hot Air Outlet**  
Do not touch or cover when hot.
- Ventilation Slots**  
Do not touch or cover.
- „AIR“ Rocker Switch (Air Pump)**  
Press the „AIR“ rocker switch, to switch on the pump for generating vapor.
- Heater Control Lamp**  
The yellow control lamp lights up during the heating process. This takes up to 5 minutes depending on the set temperature.

## FUNCTIONAL ELEMENTS

- 5 Air Filter Cap**  
Open the Air Filter Cap to replace the Air Filter underneath.
- 6 Air Filter**  
Check the Air Filter for dirt every four weeks.
- 7 „HEAT“ Rocker Switch (Heater)**  
Press the „HEAT“ rocker switch, to switch on the heater of the Vaporizer.
- 8 Temperature Regulator**  
Set the desired temperature on the setting wheel.

## FOR YOUR SAFETY

### 1 Basic Rules

- ▶ These Instructions for Use are an integral part of the Vaporizer and must be made available to the user.
- ▶ Download the current detailed version of the Instructions for Use at [www.storz-bickel.com](http://www.storz-bickel.com).
- ▶ Read these instructions carefully and completely before operating the Vaporizer.
- ▶ Always follow the instructions in these Instructions for Use when using and servicing the device.
- ▶ This Vaporizer may only be used by adults. Persons with reduced physical, sensory, or mental abilities or lack of experience and knowledge may use the Vaporizer only if they are supervised or have been instructed in the safe use of the Vaporizer and understand the resulting hazards.
- ▶ If the user is suffering from respiratory tract or lung disease, the Vaporizer may only be used after consultation and with the consent of the treating physician. Depending on their density, the vapors may irritate the respiratory tract and lungs, which may cause coughing.
- ▶ Do not operate the Vaporizer unattended.
- ▶ Store and use the Vaporizer only at a sufficient distance from heat sources (oven, stove, fireplace, etc.) and in a dry place where the ambient temperature cannot drop below +5 °C (+41 °F).
- ▶ Do not expose the Vaporizer to extreme fluctuations in ambient temperature. Condensation humidity could interfere with the functionality.
- ▶ Only use original STORZ & BICKEL accessories and spare parts.
- ▶ Do not use the Vaporizer in an inappropriate, incorrect or unreasonable manner, as the manufacturer does not accept any responsibility for any damage resulting from this.

# FOR YOUR SAFETY

## 2 Safety Recommendations

### 2.1 Risk of Injury due to electric Shock

- ▶ If the power cord is damaged, it must be replaced by a special power cord available from the manufacturer or its Customer Service.
- ▶ The Vaporizer must be switched off and disconnected from the power supply after use, before any maintenance work and during thunderstorms and in the event of faults.
- ▶ Make sure that the cable is not damaged by bending, crushing or pulling.
- ▶ Do not use multiple sockets and/or extension cables. Where absolutely necessary, use only products with a quality certificate (e.g. UL, IMQ, VDE, +S, etc.) if the stated power value exceeds the power requirements (A=Ampere) of the connected devices.
- ▶ Do not use the Vaporizer if it is damaged. Repairs may only be carried out by our Service Center. Attempting to repair damage yourself will invalidate the warranty.
- ▶ The Vaporizer must not be exposed to moisture and humidity.
- ▶ Keep the Vaporizer away from pets (e.g. rodents) and pests.
- ▶ Do not use the Vaporizer in a potentially explosive or flammable atmosphere.

### 2.2 Danger of Injury due to Burning/Vapors

- ▶ Do not touch hot surfaces of the Vaporizer (Hot Air Outlet, Filling Chamber) during operation and the cooling phase.
- ▶ Place the hot Vaporizer only on a solid surface that is heat-proof and make sure that it is not covered.
- ▶ Do not use the Vaporizer near inflammable objects such as curtains, tablecloths or paper.
- ▶ Do not close, cover or obstruct the openings of the Vaporizer during operation and the cooling phase.

## FOR YOUR SAFETY

- ▶ Before switching on the Vaporizer, ensure that all vapors, e.g. from flammable cleaning agents and disinfectants, have evaporated.

### 2.3 Risk of Injury due to various Causes

- ▶ The packaging material (plastic bags, polystyrene parts, boxes etc.) must be kept out of the reach of children as it is a potential source of danger.
- ▶ Make sure that the Valve Balloon with Mouthpiece is free of inhalable foreign bodies.
- ▶ Be careful that small parts of the Vaporizer are not swallowed.
- ▶ Be aware of the danger of strangulation through the cable.
- ▶ Be aware of the risk of strangulation and suffocation through the Balloon Tube.
- ▶ Only use the Vaporizer with the recommended herbs and plant parts to be vaporized. The use of other substances may cause poisoning or fire.
- ▶ Monitor your symptoms (e.g., cough, shortness of breath, chest pain) during or after the usage and promptly seek medical attention if you have concerns about your health.
- ▶ Use only heat-resistant, food-safe, and odor-free oven tubes or oven bags.

## OPERATION

1.



- ▶ Download current detailed version of the Instructions for Use at [www.storz-bickel.com/support](http://www.storz-bickel.com/support).
- ▶ Read the Instructions for Use before using the Vaporizer.

2.




- ▶ Grind the dried plant material in the Herb Mill.

3.



- ▶ Unscrew the Filling Chamber Cap from the Filling Chamber Housing by turning it counterclockwise.

4.

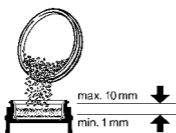

- ▶ Fill the Filling Chamber with the desired amount of shredded plant material.

5.



- ▶ Screw the Filling Chamber Cap onto the Filling Chamber Housing by turning it clockwise.

## OPERATION

6.




- ▶ Place the Filling Chamber on the hot air outlet.

7.




- ▶ Switch on the heater by pressing the "HEAT" [1] rocker switch.
- ▶ Set the desired temperature on the setting wheel [2].
- ⇒ The end of the heating process is signaled by the yellow control lamp going out.

8.

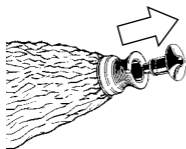

- ▶ Pull the Mouthpiece off the Valve Balloon.

9.

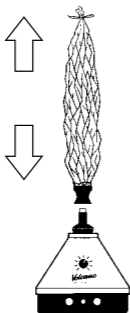

- ▶ Tighten the Valve Balloon.
- ▶ Place the Valve Balloon vertically on the Filling Chamber.



## OPERATION

10.




- ▶ Switch on the air pump by pressing the "AIR" rocker switch.
- ▶ As soon as the Valve Balloon is filled with the required quantity of vapor, switch off the air pump.

11.




- ▶ Remove the Valve Balloon from the Filling Chamber.
- ▶ Remove the Filling Chamber from the Hot Air Generator. To do this, grasp the Filling Chamber Cap by the rotary knobs.
- ▶ Place the Filling Chamber on a heat-proof surface.

12.

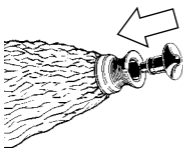

- ▶ Insert the Mouthpiece into the Valve.


13.

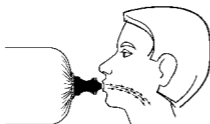

- ▶ Start the application immediately.
- ▶ Press lightly against the Mouthpiece with your lips.
- ▶ Inhale evenly for a few seconds and only half as much as possible.

## CLEANING/HYGIENE


### Explosion Drawing of the Filling Chamber Cap

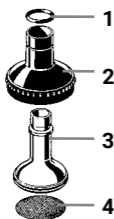

- 1 Cap Ring
- 2 Cap Housing
- 3 Cap Cylinder
- 4 Upper Screen


### Disassembling the Filling Chamber Cap

1.

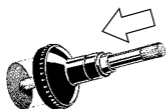

- ▶ Push out the upper Screen with the handle of the Cleaning Brush pointing downwards.


2.

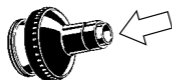

- ▶ Slide off the cap ring.
- ▶ Push the cap cylinder out downwards.


## CLEANING/HYGIENE

### Explosion Drawing of the Filling Chamber

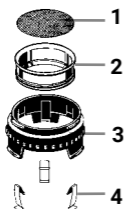

- 1 Lower Screen
- 2 Filling Chamber Cylinder
- 3 Filling Chamber Housing
- 4 3 pcs. Filling Chamber Clips

### Disassembling the Filling Chamber

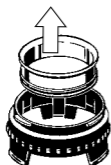

- ▶ Press the lower Screen upwards out of the Filling Chamber Cylinder.
- ▶ Press the Filling Chamber Cylinder upwards out of the Filling Chamber Housing.

## CLEANING/HYGIENE

### **Cleaning of the Filling Chamber together with the Filling Chamber Cap with rinsing water / in the dishwasher**

- ▶ Clean the Screens, the cap ring and the Filling Chamber Clips with warm rinse water, as these parts are too small for cleaning in the dishwasher.
- ▶ Clean all other parts in the dishwasher.
- ▶ Let all parts dry completely.

### **Alternative: Cleaning the Filling Chamber with Filling Chamber Cap with isopropyl alcohol**



- ▶ Soak all parts in isopropyl alcohol for a maximum of 1/2 hour.
- ▶ Rinse all parts under running water for at least one minute.
- ▶ Let all parts dry completely.

### **Assembly**

- ▶ Reassemble the Filling Chamber and the Filling Chamber Cap in reverse order.





**STORZ & BICKEL**

STORZ & BICKEL GmbH

In Grubenäcker 5-9 · 78532 Tuttlingen/Germany

